Single cell transcriptomics

Introduction to single cell RNA-seq




Why single cell RNAseq?


- Cell specific effects of treatments
- Differential gene expression between and within cell types
- Understand:
 - drug action/delivery
 - cell fate
 - cell types affected by (viral) infection
 - regulation through co-expression
 - differences in cell populations (between tissues)

• ...

scRNA-seq workflow

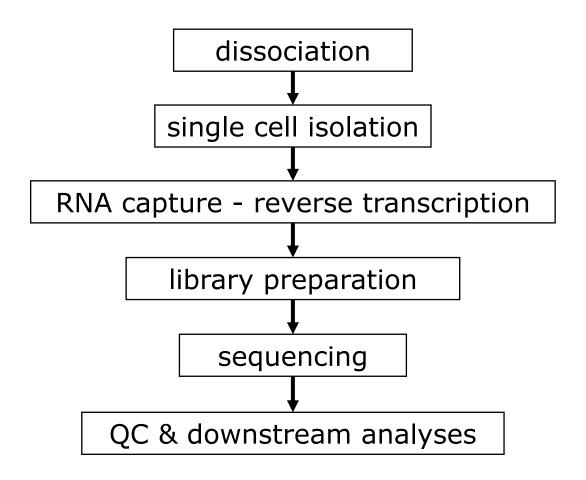

Technologies

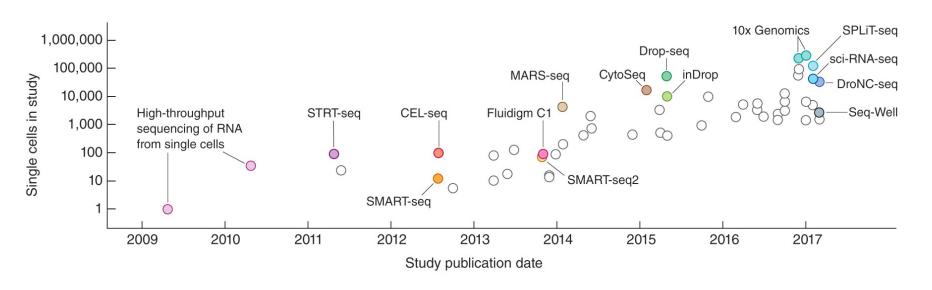

- Plate-based: SMART-seq
- Droplet-based: 10x genomics (3' kit)


Svensson V et al., Exponential scaling of single-cell RNA-seq in the past decade. Nat Protoc. 2018;13:599–604.

SMART-seq

- Use FACS for: 1 well = 1 cell
- Library preparation per cell
- Whole gene can be sequenced


6

cell lysing + reverse transcription + template switching

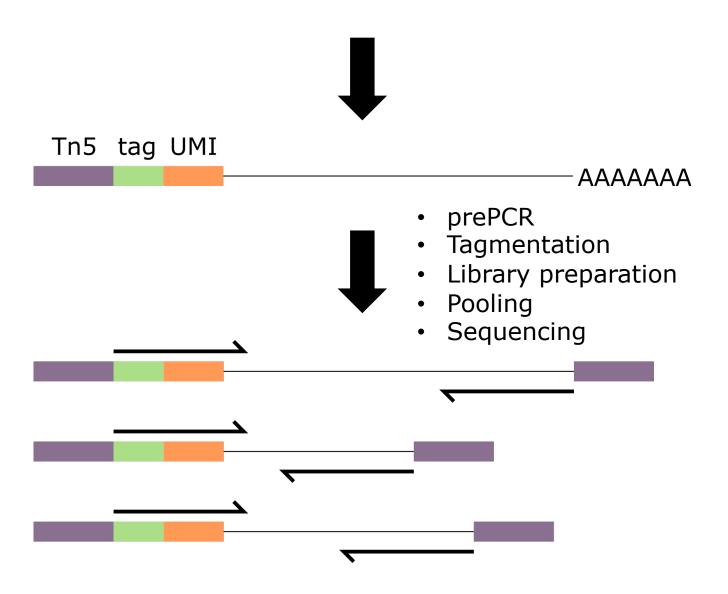



10x genomics (3' kit)


- Cell isolation within oil droplet
- RNA capture with 1 bead/cell
- Sequencing from 3' end
- Only expression

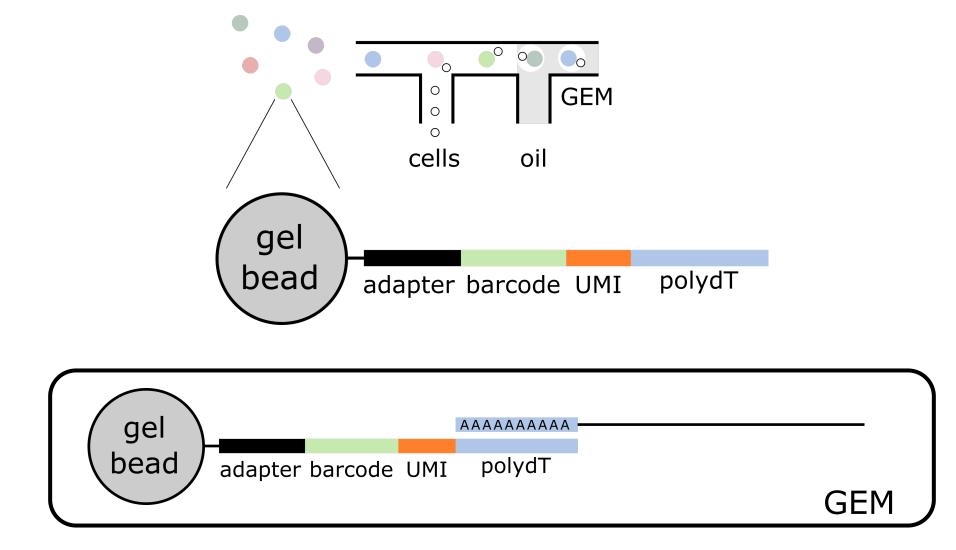

All captured **transcripts** from **single** cell: **identical** + **unique** barcode

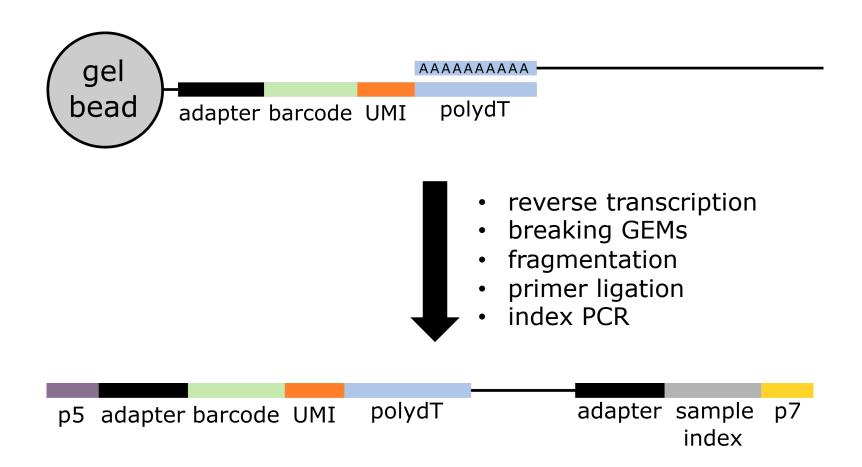

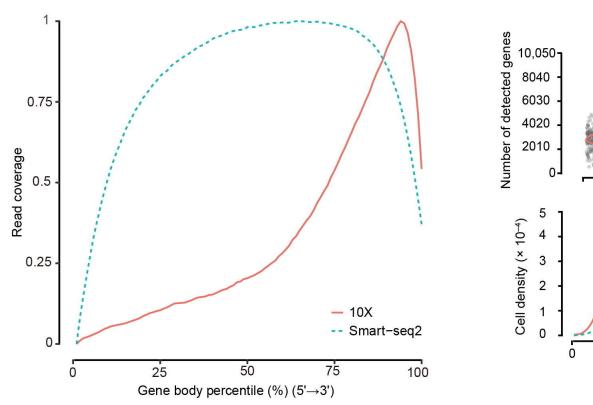

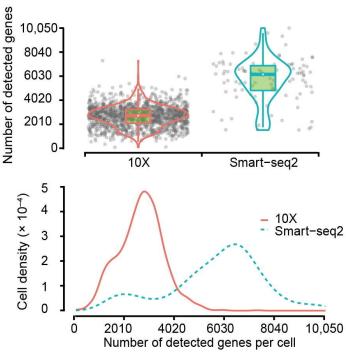

Wang X, et al. Direct Comparative Analyses of 10X Genomics Chromium and Smart-seq2. Genomics Proteomics Bioinformatics; 2021

| Droplet (10x genomics) | SMART-seq | |
|----------------------------------|---|--|
| 3' poly-A bead capture | polydT reverse transcription | |
| Strong bias 3' end | Whole transcript coverage | |
| Expression analysis | Expression + isoform analysis | |
| Low # transcripts/cell | High # transcripts/cell | |
| Investment (cell sorter) | Only FACS needed | |
| 10-100k cells | up to 1k cells | |
| 1 sample ≈ 1 library | 1 cell = 1 library | |
| Isolation by droplets - doublets | Isolation by FACS – bias to large cells | |
| Low cost/cell | High cost/cell | |

Quiz Question 2

Experimental design

- Replication, randomization and blocking
- Be aware of confounding factors, e.g.:
 - Person performing handling
 - Reagents
 - Sequencing lane/library




Record any factor for downstream correction

Further reading:


- https://doi.org/10.3389/fcell.2018.00108
- https://doi.org/10.1093/bib/bby007
- https://doi.org/10.1093/bfgp/elx035

Experimental design

Confounded design

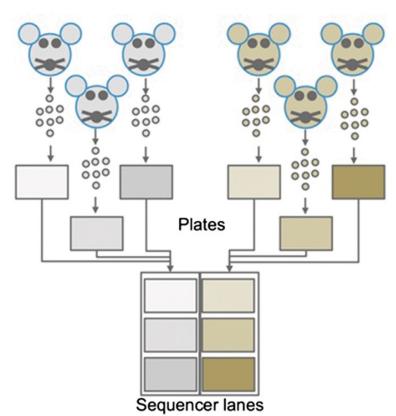

Balanced design

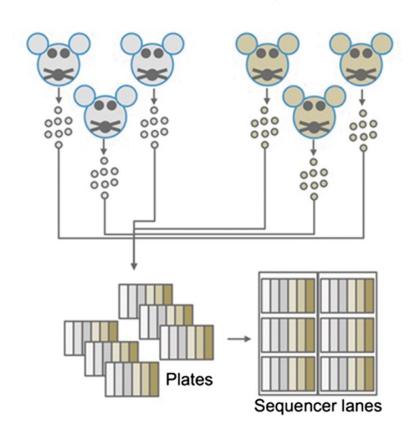


15

Baran-Gale J, et al. Experimental design for single-cell RNA sequencing. Brief Funct Genomics. 2018;17:233–9.

Question 3


Cell hashing

Solution to 1 sample = 1 library

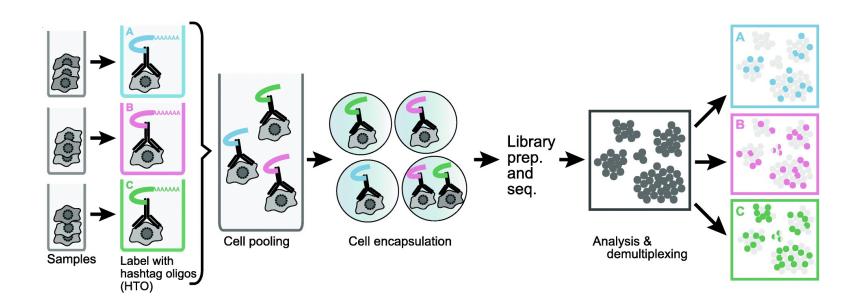

Stoeckius M et al. Cell Hashing with barcoded antibodies enables multiplexing and doublet detection for single cell genomics. Genome Biol. Genome Biology; 2018;19:1–12.

Quantify proteins

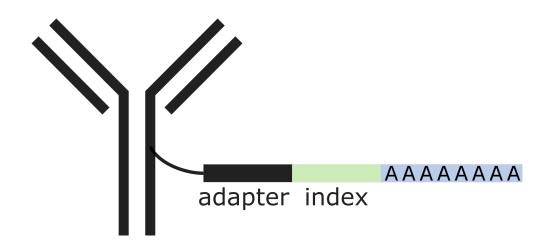

- Quantification of (cell-surface) proteins
- Together with transcriptome

Stoeckius M et al. Simultaneous epitope and transcriptome measurement in single cells. Nat Methods. 2017;14:865–8.

Question 4

single-nucleus RNA-seq


- Alternative to scRNA-seq
- For tissues difficult to dissociate
- No ribosomes -> no translation of transcription factors during processing
- Lower representation of immune cells + surface proteins

After sequencing (preprocessing)

- 1. Demultiplexing
- 2. Alignment
- 3. Quantification: # UMI/gene
- 4. Cell calling

For 10x all with cellranger




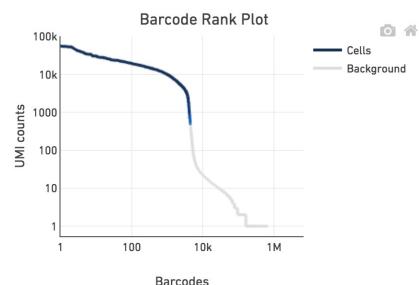

cells ->

| | ATAC-1 | CCG-1 | GCGA-1 |
|---------|--------|-------|--------|
| RPL22 | 5 | 13 | 3 |
| PARK7 | 0 | 9 | 3 |
| ENO1 | 1 | 35 | 0 |
| PLA2G2A | 0 | 0 | 0 |
| CAMK2N1 | 0 | 6 | 0 |
| CDC42 | 0 | 8 | 1 |
| C1QA | 0 | 0 | 25 |
| C1QC | 0 | 0 | 25 |
| C1QB | 0 | 0 | 29 |
| ID3 | 0 | 35 | 0 |
| RPL11 | 18 | 29 | 16 |
| CLIC4 | 0 | 4 | 0 |

Downstream analysis

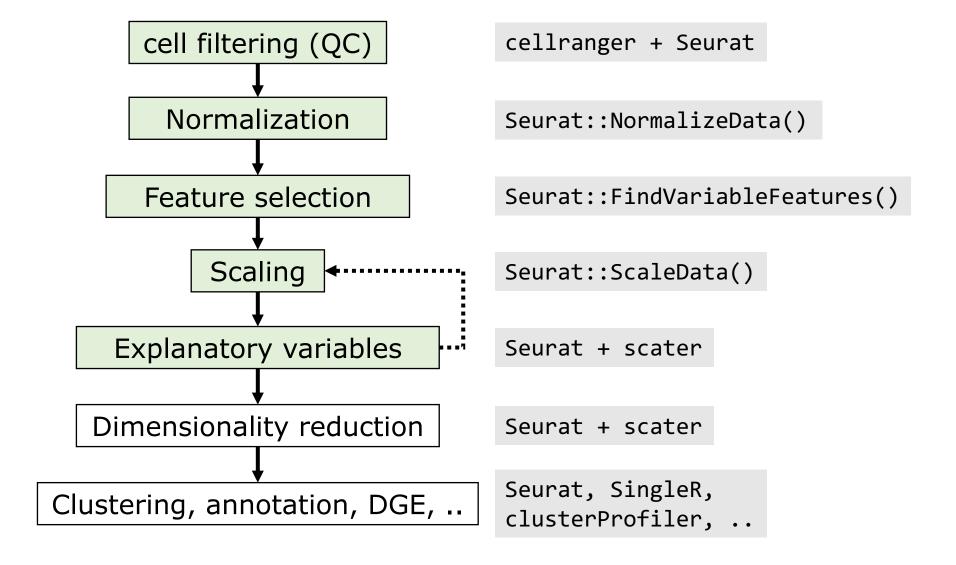

Cell filtering

- Find doublets (e.g. with DoubletFinder)
- % UMI in mitochondrial genes:
 - Points to dying/stressed cells
 - Filter by threshold (e.g. 20%)
- Number of detected genes
 - Can point to dropouts/doublets

Normalization & scaling

- Normalization (per cell): remove technical effects (i.e. library size)
- Scaling (per gene): standardize range, mean and variance

Both are mainly for the purpose of dimensionality reduction